

## Installation Instructions


Congratulations on purchasing you Aipplica Spectrometer !

Purpose of this installation Guide is to walk you through the process of connecting your spectrometer and installing necessary drivers, libraries, 3<sup>rd</sup> party or Aipplica provided, and capturing a spectra with a test source of your convenience.

You should have received a Spectrometer module along with i) an SMA905 connector 50um Core optical fiber cable with black dustcaps and ii) a USB2 (Male) to USB2 (Female) connector cable. Your Spectrometer module has a USB port for communicating with the Laptop/Desktop PC and an Optical port with SMA905 connector for Optical signal input.

Plug in one end of your USB Cable into you computer and the other end into USB2 port of the Spectrometer.

Remove the dustcap from one end of the SMA905 Optical fiber cable. Do not touch the open end with your bare hands. Connect it to the Spectrometer SMA905 port. Do not apply excessive force and only tighten it lightly. Connect the other end of SMA905 Optical cable to an optical light source or calibration source with an SMA905 port or just point it to a light source, like an LED, LED light, for coarse observation.

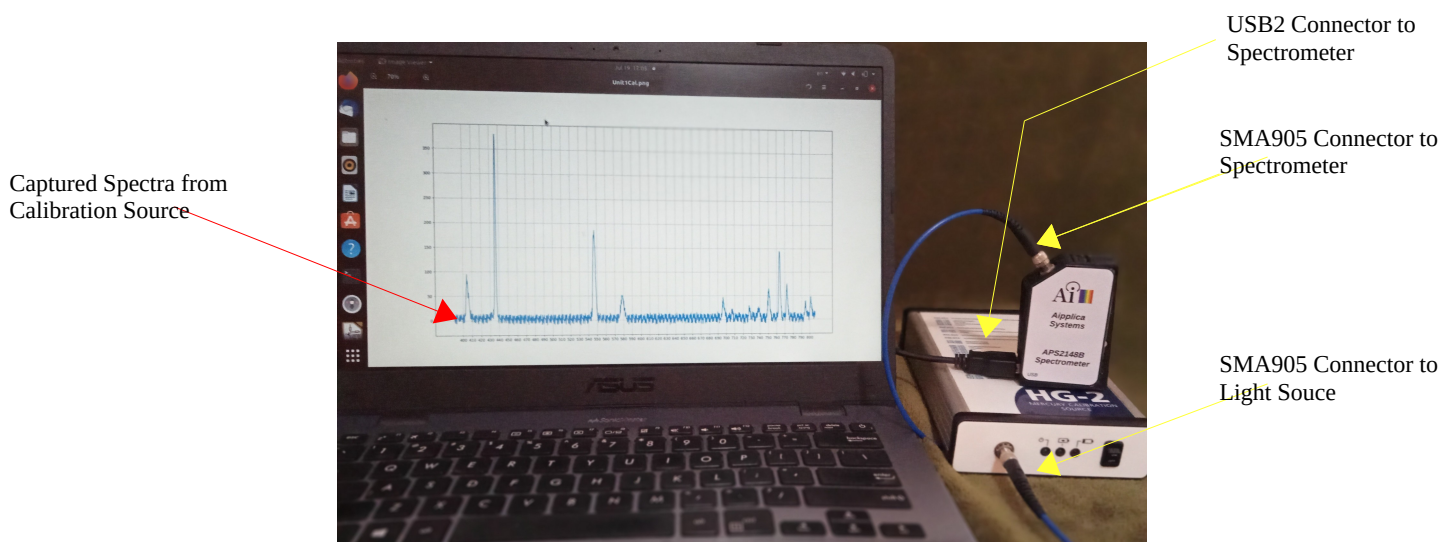

*Setup showing Aipplica Spectrometer and Calibration light Source*

From Aipplica website Download APSpec library for Python/Linux and follow the DemoA.ipynb jupyter notebook or DemoA.py to capture Spectral data.

If you are working with C language, a C library (libAPSpec.a and libAPSpec.h) with an example MyApp.c (using gnuplot) is provided. Following pages have further instructions. These examples are created with Hg-Ar calibration light source. If you have questions, you may write to [contact.aipplica.com](mailto:contact.aipplica.com) or refer to [www.aipplica.com](http://www.aipplica.com)

**NB:** Do not try to disassemble or unscrew the Unit or subject it to high mechanical stress or higher temperature than rated as it will disturb the alignment and render the calibration void and in worst case render the unit inoperative due to misalignment.

## Installation Instructions

### APSpec2148 driver python3 installation instructions.

APSpec interacts with PC through USB bus. To work with APSpectrometers in python/Ubuntu environment, you may install opensource pyftdi user space driver for popular FTDI devices. Please go through the installation section of pyftdi driver and install pyftdi package.

Alternatively you may write your own C/C++ application code and interact with Aiplica Spectrometers via FTDI's D2XX API drivers (Windows) or libftd2xx (Linux) & FTDI's MPSSE library and Aiplica Spectrometer library, libAPSpec. Go to the next section for instructions.

### For Python/Linux Environment:

#### >pip3 install pyftdi

Please note, as outlined in the pyftdi installation section, On Linux, you also need to create a *udev* configuration file to allow user-space processes to access to the FTDI devices. After you have successfully installed pyftdi.

Download APSpec2148 latest APSpec\_\*.whl from [www.aiplica.com/downloads.html](http://www.aiplica.com/downloads.html)

#### >pip3 install [APSpec-0.2-py3-none-any.whl](#)

Now you are ready to use the APSpectrometer. For plotting and viewing captured spectral data, these scripts use matplotlib. Please ensure matplotlib is installed.

#### \$pip3 install matplotlib

Invoke Python3 and follow the sequence of commands as provided in **DemoA.py** or if you use jupyter notebook, invoke jupyter notebook with **DemoA.ipynb** These files are provided in Download section on Aiplica website.

For DemoA.py executable file you may invoke the python3 executable as follows,

#### >./DemoA.py

but first read through the comments. Do block the optical input with the dustcap when collecting the dark signal and illuminate/connect the slit/fiber when collecting the light spectra at the right steps. The script will open the plot and display the captured signal. It will proceed to the next step after you close the graphics plot. Once you are familiar with the script, you may customize/modify as needed.

For DemoA.ipynb you may execute in jupyter notebook, Please read through the comments to understand the flow.

#### >jupyter notebook DemoA.ipynb

## Installation Instructions

By this point you have instantiated a spectrometer module and captured the dark signal data. You may turn on the light source of your choice (Colour LEDs, LED torch light, light bulb and capture its spectra a plot it and process the captured spectrum as you please.

### For C/C++ & Linux Environment

#### **FTDI driver installation instructions.**

APSpec uses FTDI's USB interface IC with Multi Protocol Synchronous Serial Engine (MPSSE) hardware. The MPSSE is a hardware block found in several FTDI chips which communicates with a PC over the USB interface. Applications on a PC communicate with the MPSSE in these chips using the D2XX USB drivers and the LibMPSSE library and the libAPSpec wrapper library.

Please install FTDI's libftd2xx drivers in Linux ( D2XX API drivers for Windows) Please refer to FTDI's website for additional details. The following Application Note on FTDI's website may be helpful: *AN\_220\_FTDI\_Drivers\_Installation\_Guide\_for\_Linux-1.pdf*

Next Download the LibMPSSE.zip library from FTDI's website: You will need the libMPSSE.h and libMPSSE.a files from here.

<https://ftdichip.com/software-examples/mpsse-projects/libmpsse-spi-examples/>

Download "APSpec C Library Programmer's Guide" and libAPSpec.a & APSpectra.h from Aiplica website Download sections and compile with the libraries as follows

```
>gcc MyApp.c -o MyApp.o ./libAPSpec.a ./libMPSSE.a ./libftd2xx.a -lpthread -ldl -lrt -Wall -Wextra
```

OR

```
>g++ MyApp.cpp -o MyApp.o ./libAPSpec.a ./libMPSSE.a ./libftd2xx.a -lpthread -ldl -lrt -Wall -Wextra
```

(Temporarily remove ftdi\_sio & usbserial drivers as they are not compatible with D2XX. You can restore those later with [sudo modprobe ftdi\_sio; sudo modprobe usbserial] commands)

```
>sudo rmmod ftdi_sio; sudo rmmod usbserial
```

Plug in the APSpectrometer module and Run MyApp.o

```
>chmod +x MyApp.o
```

```
>./MyApp.o
```

If you used MyApp.c Ver 0.1 from the Aiplica Website, you may see the program throw some output messages and a gnuplot of captured spectra from the light source.

```
Initialized USB Chip
Read ID Successful 7
Verification passed 0
1 1 1 2
```

## Installation Instructions

Read ID Successful

NumofCalCoeff: 03

C[0] -0.000003

C[1] 0.282167

C[2] 409.736694

Read Cal Successful

Setting Integration Time

Reading Dark Signal

Verification passed 7

Turn ON the light source and Type a character

Reading Spectral Signal

Verification passed 0

Turn OFF the light source and Type a character

Type Ctrl+d to quit..

